

Introducing

# Silicone Foam Wound Dressing

Safe n Simple focuses on product quality, patient safety, and cost containment. Simpurity wound care dressings follow evidence based medicine, facilitates moist wound healing and improves patients outcomes while lowering costs.



**Product Details:**

Size	Item#	Box/Case	HCPC
2"x2"	SNS74502	12 Pads/16 bxs	A6209
3"x3"	SNS74503	12 Pads/16 bxs	A6209
4"x4"	SNS74504	12 Pads/16 bxs	A6209
4"x5"	SNS74520	12 Pads/16 bxs	A6210
4"x8"	SNS74532	12 Pads/16 bxs	A6210
6"x6"	SNS74536	12 Pads/16 bxs	A6210
6"x8"	SNS74548	12 Pads/16 bxs	A6210
8"x8"	SNS74564	12 Pads/16 bxs	A6211

**FEATURES:**

- Core silicone oil technology provides super gentle silicone contact layer for fragile skin
- Excellent wound exudate absorption, reduces the risk of maceration and leakage
- Latex Free
- Maintains a moist wound healing environment
- No irritation to the skin upon removal of dressing
- Minimizes trauma to the wound

**INDICATIONS:**

For the management of moderate-to-highly exuding wounds such as:


- Pressure Ulcers
- Leg Ulcers
- Non-infected diabetic foot ulcer
- Donor sites
- Post-operative wound
- Lacerations and abrasions

**DRESSING STRUCTURE:**

Simpurity Silicone Foam Wound Dressings consist of a multi-layer design construction for maximum protection

Top tan color PU film backing --- be breathable, waterproof and bacteria resistance.

PU foam absorbent layer--- soft flexible foam can absorb large amount exudates.

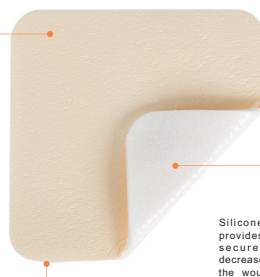

Silicone net contact layer--- provides gentle adhesion and secure fixation, especially decreases pain without damaging the wound or surrounding skin when changed..

Bottom release film